

TEVENS



Inspired by the care you deliver

Medical Device Reprocessing



Sterilization



RS120A Rotosealer

- Offers practical means for monitoring of sealing parameters
- Temperature is adjustable from 100°C to 200°C, with a tight tolerance of ±1% only, and its control is in conformance to DIN 58953-7
- Default setting for the sealing pressure is at 100N ±5%, but depends on the thickness of the material to be sealed. The upper limit is adjustable
- The display shows a warning and a sound alarms upon incorrect values

#164-RS120AV

WIPAK

ERIKING®

Visual check of the seal integrity When sealed, the film turns to darker green, so you can check that the seal is undamaged and detect any defects.

2 Superior peel of seal Pouches feature a specially designed peel-off seal

to facilitate opening. Note the correct peel-away 3 Strong performance

The seal, comprising several narrow parallel seams (Multiline), is extremely strong and durable yet peels off cleanly.

Heat-sealed upper corners reduce the risk of dust contamination when opened.

5 Easier opening Thumb notches at each end make it easier to fill and open the pack.

6 Print design safety

Indicators and text imprints are outside the actual packing area to avoid the risk of ink contamination.

Manufacture traceability

Each Steriking® package carries a lot code that states its manufacture year and month (YYMM).

Process indicators printed on the packs *EN 867-2 compliant, providing clearly contrasting* color changes

9 Exclusive materials

Made of a tough and flexible multilayer plastic film and high grammage medical grade paper.

10 Compliant with EN 868-5

Transport cartons labelled with a CE-mark. All Steriking see-through packages are manufac-tured in highly controlled Clean Room conditions to minimize the bioburden of materials.



Sterilization Packaging WPAK STERIKING®

Clear Flat Rolls and Pouches

Ideal for individual items, heavier and larger items, and smaller sets

Rolls, Flat, Clear

#164-R39 #164-R39/31	2" x 655ft 2" x 100ft
#164-R40	3" x 655ft
#164-R40/31 #164-R41	3" x 100ft 4" x 655ft
#164-R42 #164-R43	6" x 655ft
#164-R43 #164-R44	8" x 655ft 10" x 655ft
#164-R45 #164-R41/31	12" x 655ft 4" x 100ft
#164-R42/31	6" x 100ft
#164-R43/31	8" x 100ft 10" x 100ft
#164-R44/31 #164-R45/31	12" x 100ft

Pouches, Flat, Clear, Pack of 100

rack of Too	,
#164-S2 #164-S20 #164-S21 #164-S23 #164-S23	2" x 10" 8" x 16" 10" x 19½" 4" x 22½" 6" x 16" 11 ³ ⁄ ₄ " x 22½" 6" x 8" 10" x 15" 2" x 8" 3" x 12" 4" x 14" 3" x 8" 6" x 12" 11 ³ ⁄ ₄ " x 19½" 3" x 9" 3" x 9" 3" x 6" 4" x 6"

#164-S26 6" x 14" #164-S27 4¾" x 15¾" #164-S28 6" x 20½" #164-S29 10½" x 13¾" #164-S30 10½" x 17½" #164-S31 16½" x 19½" #164-S32 16¾" x 23½" #164-S33 7" x 12½" #164-S34 6⅓" x 23¾" #164-S34 6¾" x 23¾" #164-S34 6¾" x 10½" #164-S4 3" x 10½" #164-S7 8" x 10½" #164-S7 8" x 10½" #164-S7 8" x 10½" #164-S9 4" x 16"

Pouches, Self-Seal, Clear, Pack of 200

rack of 200	
#164-SS1 #164-SS10 #164-SS2 #164-SS3 #164-SS4 #164-SS5 #164-SS5 #164-SS7	3½" x 8" ½" x 8" 3.6" x 10" ½" x 22½" 5" x 10" 5½" x 15" 8" x 14" 7½" x 13" 10" x 16" 12" x 18"





Gusseted Flat Rolls and Pouches

Ideal for more voluminous items

Rolls, Gusseted, Clear

#164-RB50	3" x 1" x 330ft		
#164-RB51	4" x 2" x 330ft		
#164-RB52	6" x 2" x 330ft		
#164-RB53	8" x 2" x 330ft		
#164-RB54	10" x 2½" x 330ft		
#164-RB55	12" x 3½" x 330ft		
#164-RB56	14" x 3½" x 330ft		
#164-RB57	16" x 3½" x 330ft		

oucnes, Gussete	d, Clear
#164-B33 #164-B34	4" x 2" x 14", case of 500 6" x 2" x 16", case of 500 6" x 2" x 18", case of 500 8" x 2" x 20", case of 250 12" x 3" x 22", case of 250 4" x 2" x 12", case of 500 8" x 2" x 15", case of 250 10" x 2½" x 19", pack of 50 16⅓" x 4" x 22½", case of 250



Dupont Tyvek Flat Rolls and Pouches Ideal for low temperature sterilization

Rolls, Flat, DuPont™ Tyvek® #164-LTR40 #164-LTR42 #164-LTR42 #164-LTR43 #164-LTR43A #164-LTR44 #164-LTR45 3" x 330ft 4" x 330ft 6" x 330ft 8" x 330ft 9" x 330ft 10" x 330ft 12" x 330ft

Pouches, Flat, DuF	Pont [™] Tyvek®
#164-LTS1025	4" x 10", box 200
#164-LTS1530	6" x 12", box 200
#164-LTS1644	61/3" x 171/3", case of 600
#164-LTS1660	61/3" x 231/3", case of 600
#164-LTS2038	8" x 15", box 200
#164-LTS2538	10" x 15", box 200
#164-LTS2550	10" x 191/2", box 200
#164-LTS7520	3" x 8", box 200



Seal Integrity and Validation

Steriking® Multi Seal Test Kit

For testing the integrity of the seal made by any Rotosealer[™] units.

Kit Includes:

- Seal control sheets
- 80mL bottle of integrity test dye
- Stopwatch
- Registration card
- Instructions for use

#164-MS300

box of 250

Steriking® Seal Control Sheets

 A practical Seal Control Sheet for daily heat-sealer validation and seal quality test



• Made from the same material as the Steriking® sterilization pouches

#164-SC250

paper and laminate

box of 250

#164-TSC200

laminate and Tyvek®

box of 200

Personal Protection

Face Masks and Shields

#024-PPR-053C



#263-PG4-1201 Primagard 80® Procedure Face Masks w/ earloops, box of 50 #960-695 N-95 particulate respirator, box of 35

Full-Face Shield, clear, 13" x 8½", case of 100

#024-PPR-057A Shield, pre-assembled with frame, clear, one size, pack of 10











Bouffant Caps PRI-MED.

#263-PM4-705 SMS, fluid resistive, large, box of 250 #263-PM4-703 PSB, large, blue, box of 250



Shoe Covers

#263-PM4-110 #263-PM4-111XL



PSB, blue, bulk, case of 500 PSB, blue, X-large, dispenser box, case of 300

Gloves

■ ✓ LATEX-FREE

Nitrile Exam Gloves, Powder-Free, Chemo-Rated

#001-SC361050 X-small small medium large #001-SC361300 X-small small small medium large X-large

box of 200



■ ✓ LATEX-FREE

box of 100

■Glove Dispenser

Wire glove box holder, epoxy-coated steel.

#001-SC3003 3-box holder #001-SC3000 1-box holder



■ ✓ LATEX-FREE

box of 100

Powder-Free Stretch Synthetic Vinyl Exam Gloves

#001-SC162100 #001-SC162200 #001-SC162300 #001-SC162350 small medium large X-large



Synthetic Vinyl Exam Gloves #001-SC161100 #001-SC161200 #001-SC161300 #001-SC161350

Powder-Free



small medium large X-large

■ ✓ LATEX-FREE

Accelerator-Free Nitrile Exam Gloves,

Powder-Free, Chemo-Rated

#001-SC363050 X-small #001-SC363100 small #001-SC363200 medium #001-SC363300 large #001-SC363350 X-large

box of 100



Instruments







Commonly Replaced Instruments



Sklar—Highest quality German grade stainless steel. Designed to meet the exacting standards required for a surgical setting.



Meisterhand by Integra Miltex—German grade stainless steel. Non-glare matte finish at a competitive price.



Vantage by Integra Miltex—Japanese stainless steel. Ideal for cost-contained environments.



Operating Scissors, 51/2", Straight, Sharp/Blunt

s #444-14-1055

#162-MH5-16

V #162-V95-16

Lister Bandage Scissors, 5½"

8 #444-11-1055



Multi-Cut Utility Scissors, Serrated, 71/2"

s #444-11-1280



Dressing Forceps, 5", Serrated

s #444-19-1050

MH #162-MH6-6

1 #162-V96-6



Adson Dressing Forceps, 43/4", Serrated

1 #444-50-3147

MH #162-MH6-118

#162-V96-118

Semkin Dressing Forceps, 5", Serrated

s #444-22-4150

MH #162-MH6-104

v #162-V96-104



Adson Tissue Forceps, 43/4", 1 x 2

#444-50-3047

MH #162-MH6-120

#162-V96-120



Halsted Mosquito Forceps, 5", Straight

s #444-17-1450

#162-MH7-2

V #162-V97-2



Halsted Mosquito Forceps, 5", Curved

- **\$** #444-17-1550
- #162-MH7-4
- #162-V97-4



Schroeder Tenaculum Forceps

- s #444-91-1290 10"
- MH #162-MH30-965 93/4"
- #162-V930-965 9½"



Lucae Ear Forceps, 5½", Bayonet, Serrated

- #162-MH19-370
- #162-V919-370


Mayo-Hegar Needle Holder

- **s** #444-20-2052 51/4"
- MH #162-MH8-42 5"
- **1**#162-V98-42 5"



Baumgartner Needle Holder, 5½", Serrated

- **\$** #444-21-8005
- #162-MH8-40
- #162-V98-40

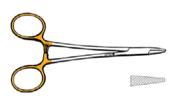

Spencer Stitch Scissors, 3½"

- **s** #444-22-2834
- #162-MH9-100
- #162-V99-100



Littauer Stitch Scissors

- **s** #444-23-1190 5½"
- ₩ #162-MH9-102 4½"
- #162-V99-102 4½"



Splinter Forceps, 41/2"

- s #444-19-3045
- #162-MH6-304
- **v** #162-V96-304



fine serrated

Bozeman Dressing Forceps, Double Curved

- s #444-80-1410 10¹/₄"
- #162-MH7-622 10½"
- **v** #162-V97-622 10½"



Foerster Sponge Forceps, 9½", Straight, Serrated Jaws

- **s** #444-87-2195
- #162-MH7-600
- **1** #162-V97-600



Mayo-Hegar Needle Holder, 6"

- **s** #444-20-2060
- #162-MH8-44
- **v** #162-V98-44



Vienna Nasal Speculum, Medium Size

- #162-20-3
- #162-MH20-3
- **v** #162-V920-3



Don't see the instrument you are looking for?

Check here for more options:

Sklar



Integra Miltex



stevens ca/gr-mdr19-2

Sterile Instruments



When it makes sense to have instruments and certain items pre-sterilized...

Scissors

#444-96-2654 Lister bandage, 5½" #444-96-2663 Littauer stitch, 3½" #444-96-2650 Iris, straight, 4½" #444-96-2651 Iris, curved, 4½" #444-96-2516 Metzenbaum, straight, 5½" #444-96-2680 Metzenbaum, curved, 7" Spencer stitch, 4½" #444-96-2596 #444-96-2522 Operating, straight, sharp/sharp, 5½" #444-96-2649G Operating, straight, sharp/blunt, 5½" Mayo, straight 5½" #444-96-2530 #444-96-2402 Knowles bandage, 5½"








#444-96-2402



#444-96-2530

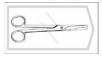



Forceps

#444-96-2708 Rochester-Ochsner, 1x2 teeth, curved, 61/4" #444-96-2532 Foerster Sponge, straight, serrated, 9½" #444-96-2554 Hartmann Mosquito, straight, 3½" #444-96-2555 Hartmann Mosquito, curved, 3½" Bozeman, S-curve, 101/4" #444-96-2620 Tissue, 1x2 teeth, 5" #444-96-2661 #444-96-2410 Fine point splinter, 4½" #444-96-2560 Kelly, straight, 5½" #444-96-2562 Kelly, curved, 5½" #444-96-2660 Dressing, serrated, 5" #444-96-2655 Adson, serrated, 43/4" Adson, 1x2 teeth, 43/4" #444-96-2656 #444-96-2626 Schroeder Tenaculum, 10" #444-96-2657 Halsted Mosquito, straight, 5" #444-96-2658 Halsted Mosquito, curved, 5" #444-96-2550 Crile, straight, 5½" #444-96-2552 Crile, curved, 5½"







#444-96-2661









#444-96-2620







#444-96-2657



#444-96-2560

#444-96-2626

Needle Holders

#444-96-2659 Webster, smooth, 51/4" #444-96-2590 Crile-Wood, 6" #444-96-2594 Mayo-Hegar, 6" Derf, serrated, 43/4" #444-96-2582 #444-96-2580 Baumgartner, serrated, 5½"







#444-96-2580

#444-96-2582

Low minimum orders!

Ask us for a quote on any items you would like pre-sterilized

DeRoyal Sterile Put-Ups

Selection includes:

- Adhesive tapes
- Bags
- Brushes
- Connectors
- Cotton-tipped applicators
- Fish hooks
- Gauze
- Instruments
- Labels
- Laboratory devices
- Lubricants
- Montgomery straps

- Plastic bottles, bowls, and cups
- Procto swabs
- Razors and blades
- Rubber bands
- Safety pins
- Sleeve protectors
- Snares
- Sponges
- Surgeon's vests
- Teflon felt and pledgets
- Tubing
- Twill tapes
- Vaginal and perineal packaging



DeRoyal

PAMA Custom Trays

Canadian manufacturer of customized sterile packs and trays.

- Fully customizable for every department in the hospital
- Ecofriendly
- Cost saving
- Fast turnaround time















Sterilization Wraps, Purple & Blue

O COI III ZU CIOI I	vviapo, i ai	
#263-PM1-2112	12" x 12" P100	case of 480
#263-PM1-2115	15" x 15" P100	case of 480
#263-PM1-2118	18" x 18" P100	case of 480
#263-PM1-2120	20" x 20" P100	case of 480
#263-PM1-2124	24" x 24" P100	case of 240
#263-PM1-2324	24" x 24" P300	case of 240
#263-PM1-2330	30" x 30" P300	case of 144
#263-PM1-2336	36" x 36" P300	case of 144
#263-PM1-2340	40" x 40" P300	case of 96
#263-PM1-2345	45" x 45" P300	case of 48







Instrument Cleaning Guide

Spaulding's Classification of Medical Equipment/Devices and Required Level of Processing/Reprocessing

Classification	Definition	Level of Processing/Reprocessing
Critical Equipment/Device	Equipment/Device that enters sterile tissues, including the vascular system	Cleaning followed by Sterilization
Semi-critical Equipment/Device	Equipment/Device that comes in contact with non-intact skin or mucous membranes but does not penetrate them	Cleaning followed by High Level Disinfectant (as a minimum requirement). Sterilization is preferred
Noncritical Equipment/Device	Equipment/Device that touches only intact skin and not mucous membranes, or does not directly touch the client/patient/resident	Cleaning followed by Low Level Disinfectant (in some cases, cleaning alone is acceptable)

Cleaning

Never keep used/dirty instruments in a dry container. Make sure BLOOD and DEBRIS don't have a chance to dry. If you are unable to process them right away use:

An Enzymatic Cleaner to help to remove organic material

#042-EBL1HC Enzyclean® instrument cleaner, pre-soak single, 1 gallon EmPower dual enzyme instrument cleaner, 1 gallon

Or a neutral pH detergent:

#042-B1HC Low suds instrument cleaner, 1 gallon

Always rinse off chemicals with distilled water.

ULTRASONIC CLEANING is usually recommended as

the best way to clean surgical instruments. #594-CSU1 Ultrasonic cleaner, 1 gallon

MANUAL CLEANING: Use stiff brushes to make sure all visible dirt is removed.

#444-10-1651-X Nylon bristles instrument brush

Use a surgical instrument cleaner

#071-002900 Organisol instrument cleaner, 4lb container, with scoop

Make sure everything looks sharp and all moving parts are working correctly.

High-Level Disinfection

Keep in mind safety, efficacy, compatibility, and the environment.

Accelerated Hydrogen Peroxide

Not recommended for instruments with brass, copper, or tungsten carbine.

Prevention—5-minute HLD for respiratory and plastics. 6-hour chemosterilant #909-11505 4L bottle

CS20—20-minute chemosterilant

#909-11405 4L bottle







Sterilization

It is recommended that all surgical instruments be autoclaved, if possible.

#144-M11-042 Midmark M11 autoclave, 11" x 18" chamber

Always use distilled water in your autoclave.

#781-AQ00334 Distilled water, 4L, case of 4 (Western Canada Only)

#320-72404 Distilled water, 4L, case of 4 (Ontario, Quebec and Eastern Canada Only)

Lubricate all moving parts.

#162-3-705 Instrument lubricant, 1 gallon

Use self-seal pouches to autoclave and store your instruments.

#164-SS1 Self-seal clear pouch, 3.5" x 8", box of 200
#164-SS2 Self-seal clear pouch, 3.5" x 10", box of 200
#164-SS4 Self-seal clear pouch, 5" x 10.5", box of 200

(Other sizes available)

Use chemical or biological indicators regularly to ensure your autoclave is working properly.

#800-1243B Steam sterilization integrator, bag of 100 #800-1262P Attest® biological indicators, box of 25 #800-116 Attest® biological indicator incubator

For instruments that can't be autoclaved or only require High Level Disinfection (HLD), use a chemosterilant soak.

Check with your Territory Manager to see what cold soak is right for you.







Environmental Infection Control

Disinfectants

#009-28+ BM-28 PLUS™ activated 2% glutaraldehyde disinfectant, 4L

#517-11-2801 MetriCide 28 high-level disinfectant, 1 gallon

Wipe down all surfaces regularly.

#907-100906722 Accel Prevention. 10" x 10" , 60 wipes anti-microbial canister, 160 wipes Isopropyl alcohol, clear, 500mL

Don't forget to clean your stethoscope and blood pressure cuffs.









Indicator Tape and Notes

#831-CCIT100 "Clean" indicator tape

100m roll

#831-CCIT100D dispenser for "Clean" indicator tape

each

#831-CCIN1000R "Not Clean" indicator notes, red

box of 4 x 250

#831-CCIN1000G "Clean" indicator notes, green

box of 4 x 250









Tip Protectors



Many other sizes and colours available!



REF Clear vented	REF Tinted Vented	REF Solid Vented	REF Solid Non-Vented	Dimension	Size	Colour*
162-3-2501C	N/A	N/A	162-3-2501	1/16" x 3/4" (1.6mm x 19mm)	1	White
162-3-2502C	N/A	162-3-2502V	162-3-2502	5/64" x 3/4" (2mm x 19mm)	2	Blue
162-3-2503C	N/A	N/A	162-3-2503	7/64" x 3/4" (2.8mm x 19mm)	3	Green
162-3-2504C	N/A	162-3-2504V	162-3-2504	1/8" x 1" (3.2mm x 25.4mm)	4	Red
162-3-2505C	N/A	162-3-2505V	162-3-2505	3/16" x 1" (4.8mm x 25.4mm)	5	Yellow
162-3-2506C	N/A	162-3-2506V	N/A	13/32 x 3/4" (10.3mm x 19mm)	6	Maroon
162-3-2507C	162-3-2507T	162-3-2507V	162-3-2507	1/16" x 3/8" x 1" (1.6mm x 9.5mm x 25.4mm)	7	Orange
162-3-2508C	N/A	162-3-2508V	162-3-2508	5/64" x 5/8" x 1" (2mm x 16mm x 25.4mm)	8	Brown
162-3-2509C	N/A	162-3-2509V	162-3-2509	1./8" x 1" x 1" (3.2mm x 25.4mm x 25.4mm)	9	Black
162-3-2510C	162-3-2510T	162-3-2510V	162-3-2510	Assorted, equal number of sizes plus 5	10	Assorted

^{*} Colour does not apply to clear vented pattern

Sterilization Carts, Totes and Bins

Pedigo Large Tables

- Stainless steel Steril-Gard construction eliminates full sub-tops where undesirable substances can collect
- 3" rubber, double ball bearing casters for silent mobility
- Sound-deadened reinforcing channels for quiet use
- Easy to clean
- Lightweight but strong

#145-SG-90-SS	w/shelf 48"W x 20"D x 34"H
#145-SG-91-SS	w/o shelf 48"W x 20"D x 34"H
#145-SG-92-SS	w/shelf 48"W x 24"D x 34"H
#145-SG-93-SS	w/o shelf 48"W x 24"D x 34"H
#145-SG-94-SS	w/shelf 60"W x 24"D x 34"H
#145-SG-96-SS	w/shelf 72"W x 24"D x 34"H
#145-SG-98-SS	w/shelf 36"W x 24"D x 34"H







Sterilization Verification and Instrument Labeling



Atttest[™] 490 Dual Auto-reader

- Intelligent display with time-to-result countdown feature
- Audible alarm sounds for improved monitoring
- Display can be viewed remotely using your smartphone or PC
- Provides clear positive (+) or negative (-) result on LCD display
- Dual Auto-reader with Any-Well technology for use with 3M[™] Attest[™] Biological Indicators 1491, 1492V, and 1295* to provide results in just 24 minutes

#800-490

#800-1322-18MM #800-1322-24MM #800-1262P #800-00135LF #800-1233LF #800-1243B #800-1250 #800-1296 #800-1298 #800-1496V #800-1491-X #800-1492V-X #242-151001EEA #175-IN-600-18 #175-IN-600-3 #175-IN-600-6 #175-IN-600-7

Comply lead-free steam indicator tape, 0.75" x 60yd roll, beige Comply lead-free steam indicator tape, 1" x 60yd roll, beige

Attest[™] biological indicators, box of 25 Comply[™] Bowie-Dick plus test, bag of 6 Comply[™] Bowie-Dick test pack, case of 5 x 6s

Comply[™] (Sterigage[™]) chemical integrator, bag of 100

Comply[™] steam chemical integrator strips, 5/8" x 8", box of 240 x 2's Attest[™] test pack rapid readout, steam, 1 cont/pack, case of 25 packs

Attest[™] rapid readout EO biological test packs, box of 25 + 25 controls Attest super rapid readout steam challenge pack, 24 controls, box of 24

Attest[™] super rapid readout indicators, gravity-displacement, box of 50 Attest™ super rapid readout Indicators, vacuum-assisted, box of 50

Instrument marking tape, white, 0.25" x 300" roll

InstroMark® instrument marking tape, black, 3/8" x 200" roll Instrument marking tape, Kelly green, 0.25" x 300" roll

Instrument marking tape, orange, 0.25" x 300" roll Instrument marking tape, royal blue, 0.25" x 300" roll

















www.stevens.ca



Atlantic Canada

Dartmouth, NS B3B 1S3

1-800-565-0765 902-468-7824 Email: ACCS@stevens.ca

Québec

101 Thornhill Drive - Unit 109 9390 boul des Sciences Anjou, QC H1J 3A9

> 1-855-660-7750 Téléc: 514-352-9214

Ontario

425 Railside Drive Brampton, ON L7A 0N8

1-800-268-0184 905-791-6143 Email: QCSAC@stevens.ca Email: ONCS@stevens.ca

Manitoba

475 De Baets Street – Unit 2 2620-61st Avenue S.E. Winnipeg, MB R2J 4M3

1-800-665-0368 1-888-640-8088

Midwestern Canada

Calgary, AB T2C 4V2

1-800-665-0368 403-640-2976

British Columbia

19051-27th Avenue. Surrey, BC V3Z 5T1

1-800-565-8444 604-585-0193 Email: MWCS@stevens.ca Email: MWCS@stevens.ca Email: BCCS@stevens.ca